wau

REN

ULTRASONIC SKIN SCRUBBER

USER GUIDE

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54/2 Fairclough Rd, Beach Haven, Auckland, Auckland,

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Table of Contents

Technical Specifications	24
Purpose	25
Device Components	25
Operating Procedure	26
Usage Method	27
Cleaning & Maintenance	30
Storage	30
Charging	30
Precautions	31
Warranty	31

Technical Specifications

Parameter	Value
Model	BPW-003
Materials	PC plastic, silicone
Dimensions (H×W×D)	150x52x24 mm
Weight	140 g
Power Supply	Lithium-polymer battery
Battery Capacity	600 mAh
Operating Voltage	3,7 V
Power	1,3 W
Ultrasonic Frequency	25 K ± 5 kHz
Charging Time	Over 2 hours
Operating Time (fully charged)	100 minutes
Smart Timer	10 minutes
Water Resistance Rating	IPX4

Purpose


The REN ultrasonic skin scrubber by WAU is the ideal device for at-home facial cleansing. It features two operating modes to provide comprehensive skincare.

In "Deep Cleansing" mode, high-frequency ultrasonic vibrations combined with positive ions help gently cleanse the skin by opening pores and removing impurities.

In "Moisturizing" mode, high-frequency ultrasonic vibrations and negative ions help soften the skin's surface, enhancing the absorption of skincare products.

The scrubber is equipped with a magnetic cap for hygienic storage and convenient transport.

Device Components

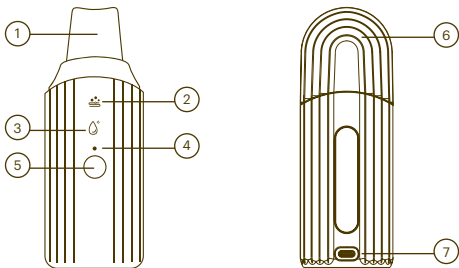1. Spatula
2. Cleansing Mode Indicator
3. Moisturizing Mode Indicator
4. Battery Level Indicator
5. Power/ Mode Button
6. Magnetic Cap
7. Charging Port

Package Contents

Device x1, USB Type-C Charging Cable x1, Storage Pouch x1
User Guide x1.


Operating Procedure

1. The device is operated using a single button. To turn it on, press and hold the Power/ Mode button. A long beep will confirm the device is powered on.
2. The Deep Cleansing mode starts automatically. The corresponding indicator will light up, and you will hear a continuous ultrasonic vibration sound.
3. To switch modes, press the Power/ Mode button again. A short beep will confirm the mode change.
4. The device will now operate in Moisturizing mode. The corresponding indicator will light up, and you will hear a pulsating ultrasonic vibration sound.
5. To turn off the device, press and hold the Power/ Mode button again. A long beep will confirm the device is powered off.
6. The device is equipped with a smart timer and will automatically shut off after 10 minutes of operation.

Usage Method

Deep Cleansing Mode

1. Before use, remove makeup and wash your face with a cleanser.
2. For best results, it is recommended to perform a cold steaming procedure using a desincrustation gel (ultrasonic cleansing gel). Apply the gel to your face, avoiding contact with the eyes. Cover with plastic film and leave for 15 minutes. After that, proceed with the ultrasonic scrubbing treatment. As an alternative to gel or lotion for ultrasonic cleansing, you may use saline solution or thermal water.
3. Press the Power/ Mode button to activate the Deep Cleansing mode.
4. Hold the device with the front side facing down (the spatula should be held at a 40-degree angle to the skin). Move the spatula along the massage lines, treating the desired areas of the face.

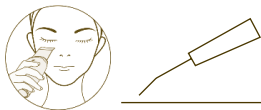5. The Deep Cleansing procedure should not exceed 5 minutes.
6. Turn off the device or switch to the next mode.
7. Rinse your face with warm water and gently clean the spatula with a dry cloth.
8. After cleansing, it is recommended to proceed to the Moisturizing mode.

Moisturizing Mode

1. After cleansing your face, evenly apply a skincare product on your choice.
2. Turn on the device and press the Power/ Mode button again to switch to Moisturizing mode.
3. Hold the device in reverse position: with the front side facing up (the spatula should lie gently on the skin rather than pressing into it). Gently glide the spatula along the massage lines, moving from bottom to top to help the product absorb fully.

4. In this mode, it is recommended to begin the procedure from the neck area, moving from the center outward.

5. Continue the massage by moving from the chin to the forehead, guiding the spatula from the center of the face toward the edges.

6. The Moisturizing procedure can last between 3 to 5 minutes. After finishing, turn off the device by pressing and holding the Power button.

7. Wipe the spatula clean and store the device properly.

Note: The recommended frequency of use for the REN scrubber depends on your skin type. For delicate or sensitive skin, cleansing is recommended no more than once a week. For oily skin, the device can be used up to twice a week.

Cleaning & Maintenance

After each use, clean the spatula and the device with a damp cloth, then wipe dry with a soft, dry cloth.

The scrubber has an IPX4 water resistance rating, providing protection against splashes. However, the device must not be submerged in water.

Storage

- Store the device in a dry place, protected from dust.
- Keep the device away from heat sources.
- Store the device out of reach of children.

Charging

The battery level is indicated by an LED light on the front of the scrubber.

When the battery is sufficiently charged, the indicator remains steadily white.

When the battery is low, the indicator blinks, and the device will automatically turn off.

To charge, connect the Type-C end of the charging cable to the device's charging port and plug the USB end into an adapter.

During charging, the LED indicator lights up red. Once fully charged (100%), the indicator turns white and stays on continuously.

Precautions

1. This device is intended for individual use only. Do not share it with others.
2. Use the device strictly as intended, adhering to the recommended duration and frequency of use.
3. Use caution when applying the device near areas with metal implants or prostheses.
4. Do not use the device in areas affected by venous or arterial disorders, or for patients with thrombophlebitis or a history of blood clotting disorders.
5. Do not drop the device.
6. Do not disassemble, modify, or repair the device.

Warranty

The REN ultrasonic skin scrubber by WAU is covered by a 1-year warranty. During the warranty period, the device is eligible for replacement or repair in the event of any technical defects. This warranty does not apply to malfunctions caused by improper use, transportation damage, or unauthorized repairs performed by the user.

REN

ULTRASONIC SKIN SCRUBBER

Model: **BPW-003**

Material: **PC plastic, silicone**

Dimensions (H×W×D): **150×52×24 mm**

Weight: **140 g**

Voltage: **3,7 V**

Power: **1,3 W**

Made in: **PRC**

Authorized manufacturer representative: PERFECT SQUARE
LIMITED 54/2 Fairclough Rd, Beach Haven, Auckland,
Auckland, 0626, New Zealand
info@wauglobal.com

Important Notice: To ensure continuous product
improvement, the manufacturer reserves the right to
modify design and technical specifications. The latest
version of this manual is available at: wauglobal.com