wau

MIGA PRO
PORTABLE WATER FLOSSER
USER GUIDE

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RoHS

Table of Contents

Technical Specifications	28
Purpose	28
Device Components	29
Operating Procedure	30
Usage Method	32
Cleaning & Maintenance	33
Storage	33
Charging	34
Malfunctions and how to fix them	34
Precautions	35
Disposal	36
Warranty	36


Technical Specifications

Parameter	Value
Model	RM-MINIO18
Rated Supply Voltage	5 V
Rated Operating Voltage	5 W
Power Source	Lithium Battery
Battery Capacity	2200 mAh
Water Resistance Rating	IPX7
Modes	Normal / Pulse / Soft
Water Tank Capacity	180 ml
Operating Time	Refer to the Operating procedure section
Время зарядки	4 часа

Purpose

The MIGA Pro portable water flosser by WAU is a portable oral care device. It features 3 operating modes, each with 3 adjustable speeds, as well as a sterilization mode. The device includes an enlarged water reservoir and comes with 3 different nozzles, ensuring effective removal of food debris, dental plaque, and thorough cleaning of interdental spaces and periodontal pockets. The device has a foldable design, making it convenient to carry while traveling or to use for the care of orthodontic appliances.

Device Components

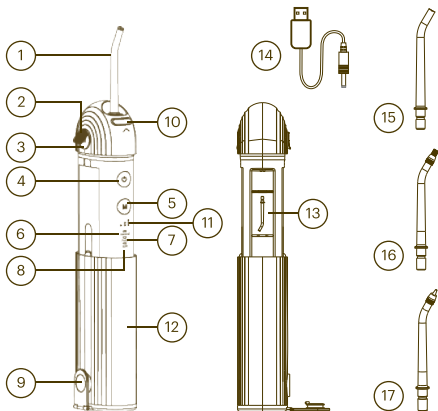1. Nozzle

2. Charging Port Cover

3. DC Charging Port

4. Power Button

5. Mode Selection Button

6. Standard Mode

7. Pulse Mode

8. Soft Mode

9. Water Tank Cover

10. Nozzle Release Button

11. Water Pressure Indicator

12. Water Tank

13. Nozzle Storage Compartment

14. DC Charging Cable


15. Standard Nozzle

16. Orthodontic Nozzle

17. Periodontal Nozzle

Operating Procedure

1. Attach the nozzle until it clicks into place.
2. Open the reservoir and fill it with water.
3. Close the tank lid tightly to prevent it from opening.
4. Briefly press the power button to turn on the device.
5. Briefly press the power button again to switch cleaning modes.
6. Press and hold the mode button for 3 seconds to change the speed.
7. Select the desired mode and speed, place the irrigator in your mouth, and press and hold the power button for 3 seconds to start cleaning.
8. After use, turn off the device by pressing and holding the power button for 3 seconds. The device will remember the last used mode and automatically start in that mode the next time it is turned on.
9. Press and hold the nozzle release button to remove the nozzle carefully. Open the nozzle storage compartment and place the nozzle inside.
10. Remember to clean regularly: detach the reservoir and wash it thoroughly.

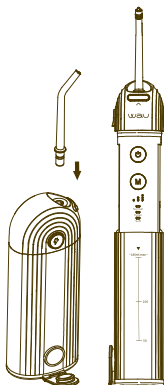[Note]: Use the gentle mode during the initial period of use. Mild bleeding during the first few cleanings is normal and should disappear within one week. If bleeding persists for more than two weeks, discontinue use and consult a dentist.

UV light sterilization function

Automatic function: when water tank is closed to main body, the UV light will automatically work for 10 minutes and then stop. Three mode indicators will flash to indicate when the UV light is sterilizing.

Reset function

Within 10 seconds of input charging, press power button 5 times to restore factory settings (all strength lights flash 5 times 1HZ), that is, restore to standard 1, weak strength. Press any button for normal operation.

Pulse Frequency on each mode and speed

Mode	Speed	Pulse Frequency
Normal	Weak	1440 / min.
	Medium	1920 / min.
	Strong	2400 / min.
Pulse	Weak	1440 / min. $\pm 10\%$
	Medium	1920 / min. $\pm 10\%$
	Strong	2400 / min. $\pm 10\%$
Soft	Weak	960 / min. $\pm 10\%$
	Medium	1440 / min. $\pm 10\%$
	Strong	1920 / min. $\pm 10\%$

Operating Time on each mode and speed

Mode	Speed	Pulse Frequency
Normal	Weak	196 min.
	Medium	72 min.
	Strong	59 min.
Pulse	Weak	133 min.
	Medium	103 min.
	Strong	84 min.
Soft	Weak	140 min.
	Medium	96 min.
	Strong	72 min.

Usage Method

a. Use the standard multifunctional nozzle. Press the power button, select the desired mode and speed. Then position the nozzle so that the water nozzle is perpendicular to the area being cleaned, aiming at the teeth. [Note]: Sensitive gums tend to bleed easily. Do not direct the jet directly at the gums!

b. Begin cleaning with the back teeth and move towards the front teeth. Continue this way along the entire dental arch, pausing briefly at each tooth until the gum line and interdental spaces are cleaned. If you started from the outer side of the teeth, repeat the procedure for the inner side. This routine should be performed for both upper and lower teeth. [Note]: For people with sensitive teeth, it is recommended to use the Gentle mode.

c. To clean periodontal pockets, use the periodontal nozzle. Press the power button, select the appropriate mode and speed. Position the nozzle at an angle so that the tip is directed at the gum line and into the periodontal pocket. Perform cleaning following the instructions described above.

d. The orthodontic nozzle is designed for cleaning braces. In addition to the water jet outlet, it features a brush with stiff bristles at the tip, which effectively cleans the brackets by removing plaque and food debris.

Cleaning

1. After use, turn off the irrigator and remove the nozzle.
2. Thoroughly wash, dry, and store it in a dry place.
3. This device has an IPX7 waterproof rating and can be rinsed directly under running water.

Storage

a. Nozzle

Press the button to release and remove the nozzle. Rinse it with water, then dry thoroughly with a cloth and store it in the nozzle storage compartment. It is recommended to replace the nozzle every 6 months. If the tip of the nozzle becomes deformed, replace it immediately.

b. Device Body

Wipe off dirt with a cloth dampened with mild detergent, then dry with a clean cloth. If the body is heavily soiled, rinse it with water. Do not use solvents, gasoline, alcohol, or similar substances. Do not wash in water hotter than 40 °C. After cleaning, be sure to dry the device body thoroughly.

c. Water Filling Compartment

Clean regularly to prevent leaks.

d. Removable Water Reservoir

Open and rinse the removable water reservoir regularly.

Clean the inside of the reservoir with a brush or similar tool, then wipe it dry with a cloth. If the device will not be used for an extended period, do not leave water in the reservoir.

Charging

1. Low electric quantity reminder: when the battery is lower than 10% during work and power-on, the LED (3 modes) light will flash to remind that the product needs to be charged.
2. Turn off the device and retract the water reservoir.
3. Open the charging port cover by gently pressing the lower part of the cover.
4. Insert the DC charging cable connector into the charging port.
5. Plug the adapter into a power outlet and make sure charging has started — during charging, the mode indicators will flash in sequence.
6. A full charge takes 4 hours. Operating time on a full battery may vary depending on the selected mode and speed (see the Operating Instructions section). When the device is fully charged, all three-mode indicators will remain steadily lit.

Malfunctions and how to fix them

Malfunction	Reason	How to fix
Not working	A new device was purchased or it hasn't been used for more than 3 months	Charge continuously for more than 4 hours
Low water pressure	a. Nozzle is clogged with dirt b. Nozzle is deformed	Replace the nozzle
Works only for a few minutes after charging	a. Insufficient charging time, low battery level b. Battery life has expired	a. Charge continuously for more than 4 hours b. Contact customer service
Battery indicator does not light up during charging	a. Charging cable is not inserted correctly b. Power circuit is damaged	a. Check that the cable is connected properly b. Contact customer service

Precautions

1. Keep the charging port dry during charging.
2. Do not connect the device to a power source with a voltage different from the specified one.
3. Do not touch the charger with wet hands.
4. Do not use the device if it is damaged or malfunctioning.
5. Do not aim the water jet directly at the tongue, ears, or other parts of the body — water pressure may cause injury.
6. If you have dental conditions, use this device under the supervision of a dentist.
7. This device can be used together with mouthwash or other products recommended by your dentist.
8. Do not disassemble the device or insert foreign objects into any openings unless you are a qualified technician.
9. Do not operate the device in harsh outdoor environments.
10. Do not fill the water tank with iodine, bleach, or water-insoluble aromatic oils, as these substances may damage the device and shorten its lifespan.
11. Remove all oral jewelry before using the device.
12. Do not use the device if there are visible wounds on the tongue or inside the mouth.
13. If prescribed antibiotics for treatment, do not use this device unless specifically advised by your doctor.
14. Do not direct the water jet into the ears to avoid injury.
15. Keep the device out of reach of children under 3 years old.
16. Store the device in a dry, clean place away from direct sunlight and heat sources.
17. Before storing the device for a long period (over 3 months), fully charge the battery.
18. Do not operate the device without water in the tank, as this may cause malfunction.

[Notes]:

1. Do not use hot water above 40 °C to clean the water tank. To prolong the life of the tank, avoid deforming it due to exposure to high temperatures. It is recommended to use room temperature water below 40 °C or running tap water for cleaning this device.

2. Do not add toothpaste, salt, hard mouthwash tablets, or other chemicals to the water to avoid clogging and damaging the device.
3. Before adding mouthwash or fresheners, check that they do not contain isopropyl myristate, which can damage the device and affect its normal operation.
4. Do not use thinners, alcohol, or other harsh cleaning agents, as this may cause discoloration and malfunction of the device.

Disposal

Do not dispose of the device with household waste. Instead, take it to a specialized recycling center to protect the environment.

[Battery Handling]

The built-in battery contains substances that may pollute the environment. When disposing of the irrigator, be sure to remove the battery and recycle it at an official recycling center.

Do not disassemble the battery; only remove it prior to disposal. Before removing the battery, ensure it is fully discharged.

[Battery Removal]

1. Disconnect the irrigator from the power source.
2. Turn off the irrigator.
3. Open the housing and remove all parts.
4. Cut off the positive and negative battery terminals.
5. Remove the battery.

Warranty

The MIGA Pro Portable water flosser by WAU comes with a 1-year warranty from the date of purchase. If technical defects are found within the warranty period, the device will be replaced or repaired.

This warranty does not cover malfunctions caused by improper use or transportation, or devices that have undergone self-repair.

MIGA PRO

PORTABLE WATER FLOSSER

Model: **RM-MINI018**

Materials: **ABS plastic, PC plastic**

Dimensions (HxWxD): **172 x 63,5 x 41,7 mm**

Weight: **260 g**

Voltage: **5 V**

Power: **5 W**

Made in: **PRC**

Authorized manufacturer representative: PERFECT SQUARE LIMITED 54/2
Fairclough Rd, Beach Haven, Auckland, Auckland, 0626, New Zealand
info@wauglobal.com

Important Notice: to ensure continuous product improvement, the manufacturer reserves the right to modify design and technical specifications. The latest version of this manual is available at: wauglobal.com

RoHS