wau

LED SET

PHOTOTHERAPY FACE & NECK MASK

USER GUIDE

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Table of Contents

Technical Specifications	17
Purpose	18
Device Components	18
Operation Procedure	20
Usage Method	20
Precautions	23
Cleaning & Maintenance	24
Storage	24
Charging	25
Disposal	26
Warranty	26
Frequently Asked Questions	26

Technical Specifications


Parameter	Value
Model	K9C-1KB/ FM-33
Materials	ABS Plastic, Silicone
Dimensions (HxWxD)	Face Mask: 190x679x5 mm, Neck Mask: 360x213x7.3 mm
Weight	Face Mask: 380 g, Neck Mask: 210 g
Power Source	Lithium-ion Battery
Battery Capacity	1500 mAh
Operating Voltage	3.7 V
Power	8 W
Charging Time	4 hours
Operating Time	30 minutes (when both masks operate simultaneously at maximum light intensity)
Smart Timer	15 minutes
Modes	Red/ Blue/ Green/ Yellow/ White Light

Purpose

The LED SET by WAU offers a comprehensive phototherapy solution with two masks designed to treat three areas: face, neck, and décolletage. Both masks connect to a single remote control, allowing you to perform light therapy sessions on multiple areas within one treatment.

Utilizing advanced phototherapy technology, these masks emit light waves of varying lengths to stimulate restorative processes in the skin. This helps to smooth wrinkles, combat pigmentation, acne, couperose, and rosacea. Regular use of the LED masks will help maintain skin tone and fight against signs of aging and external skin imperfections.


Device Components

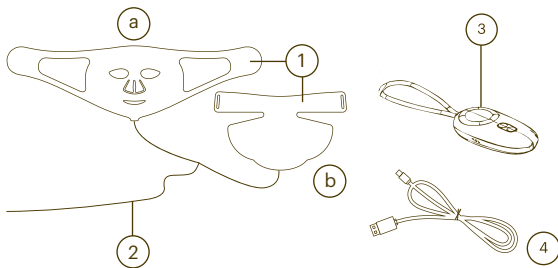1. Phototherapy Masks:
 - a. Face Mask
 - b. Neck and Décolleté Mask
2. Type-C Cable (dual output)
for connection to the remote control
3. Remote Control
4. USB Type-C Charging Cable

Remote Control Overview

1. Power On/Off Button
2. Pause Button
3. Mode Selection Buttons
4. Light Intensity Adjustment Buttons
5. Timer Display
6. Charging Indicator
7. Mode Indicator (Color)
8. Light Intensity Level Indicator

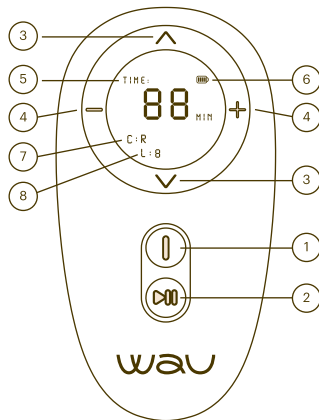Package Components

- Face Phototherapy Mask x1
- Neck & Décolleté Phototherapy Mask x1
- Remote Control x1
- Type-C Cable for Connecting the Remote to the Masks x1
- USB Type-C Charging Cable x1
- Storage bag x1
- User Manual x1

Operating Procedure

1. Connect the masks to the remote control using the Type-C cable. Insert each cable end into the mask ports until you hear a distinct click.

You may use both masks simultaneously or connect and use them one at a time.

2. Press and hold the power button on the remote control for 2 seconds to turn on the devices.
The masks will automatically start in red light mode at intensity level 1 (display shows: C: R, L: 1).

3. Briefly press (1 second) the mode selection buttons (^/v) to switch between:

Blue light mode (C: B)

Green light mode (C: G)

Yellow light mode (C: RG)

White light mode (C: RGB)

4. Briefly press (1 second) the light intensity buttons (-/+) to adjust the intensity level from 1 to 3 and back (display shows: L: 1 / 2 / 3).

5. Briefly press (1 second) the pause button to pause the operation; press it again to resume.

6. To power off, press and hold the power button for 2 seconds.

Usage Method

[Important] Perform a sensitivity test before first use:

- Place the mask on the inner side of your arm.
- Turn it on and run a full 15-minute phototherapy session.

- If you experience discomfort or a burning sensation, stop using the mask immediately.
- Wait around 6 hours. If redness or itching appears on your skin, the LED SET is not suitable for your use.

1. Connect the remote control to the masks.
2. Put on the face mask and the neck & décolleté mask.
3. Turn on the masks, select the desired light mode and intensity level.
4. Enjoy your phototherapy session.

Each mode runs for 15 minutes, after which the device will turn off automatically.

Switching modes resets the timer.

[Important]: Ensure the total duration of use in one session does not exceed 15 minutes, even if switching between modes.

5. After your session, apply your skincare product (serum or cream), as the skin's absorption ability increases after phototherapy.
6. Use the masks regularly 3–4 times per week for at least 4 weeks to achieve visible results.
7. After achieving the desired effect, reduce use to 1–2 sessions per week for maintenance.

Mode Features

Red Light (660 nm)

Red light has the longest wavelength, allowing it to penetrate the skin more deeply than other visible colors. It acts on fibroblasts, stimulating collagen production and promoting skin regeneration. Red light therapy accelerates wound healing, tightens the skin, and reduces visible signs of aging such as wrinkles.

Blue Light (460 nm)

Blue light has the shortest wavelength among visible colors. It acts at the epidermal level and has antimicrobial properties. It is effective in treating acne by destroying the bacteria that cause breakouts.


Green Light (525 nm)

Green light is absorbed by the epidermis and dermis, triggering a cellular response that leads to skin repair and healing. It helps reduce redness, soothes sensitive skin, and combats couperose by minimizing the appearance of capillaries and supporting healthy blood flow.

Yellow Light (Red + Green = Yellow)

Yellow light in this device is created by combining red and green light spectra. This polychromatic light therapy produces more comprehensive results: it stimulates blood circulation, brightens skin tone, and reduces pigmentation.

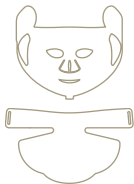White Light (Red + Blue + Green = White)

The mixed wavelength of white light reaches the dermis and boosts cellular metabolism, enhancing the skin's absorption capacity. After using the white light mode, your skin can better absorb active ingredients in skincare products. For optimal results, apply your skincare product immediately after the session.

Precautions

- LED phototherapy masks are intended for home use only. Use them strictly according to the methods described in this manual. Do not use the device for other purposes, as this may result in product damage or personal harm.
- If any part of the mask is damaged, discontinue use immediately and contact an authorized service provider. Any self-repair, disassembly, or misuse voids warranty obligations.
- Do not use the masks while charging. Avoid touching the mask, charging cable, or USB adapter with wet hands during charging.
- Stop use immediately if you feel discomfort or experience any side effects (such as skin redness, itching, or allergic reactions).
- Not suitable for infants, individuals with epilepsy, thyroid disorders, or those taking thyroid medication.
- Do not use on damaged, injured skin or post-surgical wounds.

Cleaning & Maintenance

1. Disconnect the device from any power source before cleaning.
2. Do not spray water directly onto the device.
3. Wipe the surface of each mask twice, for three minutes each time, using a 70–80% ethanol solution. Then wipe with a damp or dry cloth.
4. Do not use flammable substances such as benzene, solvents, butyl acetate, etc., for cleaning.
5. Do not use other cleaning agents, as they may damage the masks.

Storage

1. Clean the surface of the masks after each use.
2. Store the device in a dry room at room temperature.
3. Avoid storing the device in dusty areas, direct sunlight, or humid environment.
4. The masks do not contain any user-serviceable parts.

Charging

The battery level of the LED SET is displayed on the remote control.

-  At 100% charge, all four indicators in the battery icon are lit.
-  From 100% to 75%, three indicators remain steadily lit.
-  From 75% to 50%, two indicators are lit.
-  From 50% to 25%, one indicator remains steadily lit.
-  From 25% to 0%, one indicator blinks continuously.

To start charging, connect the Type-C end of the charging cable to the remote control and the end with the USB connector to the adapter.

-  At 0–25%, one indicator blinks.
-  At 25–50%, one indicator is solid, and the second blinks.
-  At 50–75%, two indicators are solid, and the third blinks.
-  At 75–100%, three indicators are solid, and the fourth blinks.
-  At 100%, all four indicators stay lit as long as the remote is connected to a power source. Once disconnected, the display turns off automatically.

- If the supply voltage drops below 3.2 V, the charging indicator blinks three times, and the remote control automatically powers off.
- Please use only the original USB cable provided in the kit for charging. Do not connect the mask directly to the USB charging cable.
- Using an adapter with less than 5V output will not charge the device. Using an adapter with more than 5V output may damage the product.
- This device includes over-voltage protection for voltages above 6 V and below 20 V.
- Do not use the product while charging.

Disposal

1. Before disposing of the device, the battery must be removed.
2. The device should be powered off before removing the battery.
3. After removal, dispose of the battery in accordance with local regulations.

Warranty

The WAU LED SET Facial and Neck Phototherapy Masks come with a 1-year warranty. If any technical defects are discovered during the warranty period, the product is eligible for repair or replacement. This warranty does not cover damage caused by misuse, improper transportation, or unauthorized repairs.

Frequently Asked Questions

Q: Are these masks waterproof?

A: No, the masks are not waterproof. They should only be wiped with water or ethanol. Do not wash them.

Q: My skin feels dry after using the masks. What should I do?

A: Light wavelengths penetrate into the dermis and stimulate cellular metabolism, which increases the skin's absorption capacity. As a result, your skin may need more hydration and nourishment. After each session, apply sufficient moisturizers or serums.

Q: Do the masks heat up during use?

A: No noticeable heating occurs when using red, blue, or green light. However, in mixed light modes (yellow/ white), the temperature of the mask may rise up to $38 \pm 1^{\circ}\text{C}$.

Q: Will I see results after the first use?

A: Visible effects are minimal after 1–2 sessions, since the skin's natural cycle is approximately 28 days. For noticeable results, use the masks consistently for at least 1 month.

Q: What's the recommended duration for one LED SET session?

A: Each session should last no longer than 15 minutes. Exceeding this may harm the skin. Always monitor your session time.

Q: Can the masks be bent or folded?

A: No. Bending or folding the masks may damage the internal circuitry.

LED SET

FACIAL AND NECK PHOTOTHERAPY MASKS

Model: **K9C-1KB/ FM-33**

Material: **ABS plastic, silicone**

Dimensions (HxWxD):

Face mask – 190 x 679 x 5 mm

Neck mask – 360 x 213 x 7.3 mm

Weight:

Face mask – 380 g

Neck mask – 210 g

Voltage: **3.7 V**

Power: **8 W**

Made in: **PRC**

Authorized manufacturer representative: I.P. Umantsev A. N. Leningradsky Prospekt 37, Office 14–14, Moscow, 125167, Russia

info@wauglobal.com

Important Notice: to ensure continuous product improvement, the manufacturer reserves the right to modify design and technical specifications. The latest version of this manual is available at: wauglobal.com