wau

MORI 2.0

SMART SONIC TOOTHBRUSH

USER GUIDE

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Table of contents

Technical Specifications	21
Purpose	22
Device Components	22
Operation Procedure	23
Usage Method	27
Cleaning & Maintenance	29
Storage	33
Charging	33
Precautions.....	35
Disposal	35
Warranty.....	35

Technical Specifications


Parameter	Value
Model	RM-T40
Power	5 W
Operating Voltage	3,7 V
Power Source	Lithium-ion Battery
Battery Capacity	3000 mAh
Water Resistance Rating	IPX7
Modes & Oscillation Frequency	Sensitive (28 000 osc/min) Clean (37 000 osc/min) White (41 000 osc/min) Polish (36 000 osc/min) Massage (35 000 osc/min)
Smart Timer	2/ 3/ 4 minutes / unlimited
Operating Time	Sensitive Mode — 1200 min Clean Mode — 250 min White Mode — 400 min Polish Mode — 500 min Massage Mode — 400 min
Charging Time	8 hours

Purpose

The WAU MORI 2.0 Electric Sonic Toothbrush effectively and gently cleans teeth while removing plaque and tartar. It's operation is based on high-frequency vibration technology: cleaning occurs through rapid bristle oscillation activated by sonic waves (high-frequency vibration).

The MORI 2.0 features 5 intensity levels per mode, a memory function, a quadrant timer, 3D oral mapping to evaluate cleaning effectiveness, a pressure-sensing indicator, brush head replacement reminders, and a UV sanitization case — all designed to ensure high-quality and comfortable oral hygiene.

Device Components

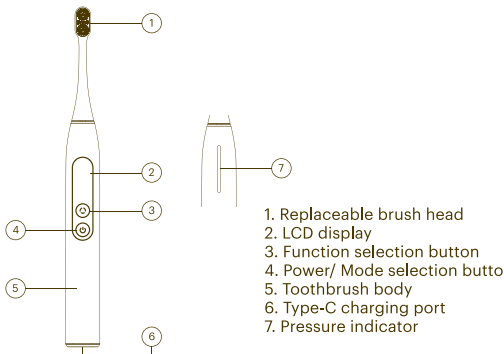Package Contents

Toothbrush body ×1, Brush heads: Sensitive ×1, Whitening ×1, Cleaning ×1, Type-C charging cable ×1, Storage pouch ×1, Instruction manual ×1, UV sanitization case ×1.

Operation Procedure


1. The MORI 2.0 toothbrush features palm-activated auto-start functionality: it powers on automatically when you grip the body. **Important:** during initial use, auto-start is disabled to prevent accidental activation during transit. To enable auto-start and 3D oral mapping, you must first activate the gyroscope by following these instructions.
2. Briefly press the power/mode selection button to turn on the brush.
3. Navigate menus using the function selection button, confirm selections with a short press of the power/ mode button.

4. Scroll to the "Settings" section and enter it. Within the menu, select "3D OFF". Activate 3D mapping by pressing the power button once.

5. The MORI 2.0 supports two language options: Russian and English. To change the language, navigate to "Settings" select the "Language" section, and briefly press the function selection button to set your preferred language — the toothbrush is now ready for further configuration and use.

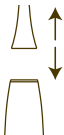
6. The kit includes 3 types of brush heads:

- For sensitive teeth, with soft bristles;
- For whitening, with firm bristles;
- For cleaning, with medium-firm bristles.

Select the required head, align it directly above the motor shaft tip (transmission shaft), and install it.

Important: Do not twist the head during installation. Note that a small gap remains between the toothbrush body and head is normal.

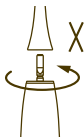MORI 2.0 brush heads feature NFC sensor chips: during initial installation, the brush automatically selects the corresponding cleaning mode — the sensitive head activates Sensitie Mode; the whitening head activates White Mode; the cleaning head activates Clean Mode.

If the display shows an "Error" symbol with a crossed-out brush icon, this indicates an incompatible brush head has been installed. In this case, the toothbrush cannot be activated.

6. After installing the brush head, continue configuring your MORI 2.0 toothbrush as instructed. Navigate to the "Settings" section.

5. Select the "Brushing Time" section. Set the duration: 2, 3, or 4 minutes, or "Unlimited" by pressing the function selection button. The MORI 2.0 is now ready for further configuration and use.

7. Return from the "Settings" section to the main menu by selecting the "Back" option.

8. Select your cleaning mode in the menu:

Sensitive Mode is recommended for first-time users of electric sonic toothbrushes and those with increased tooth sensitivity.

Clean Mode is designed for daily oral hygiene routines.

White Mode effectively removes stains and addresses tooth discoloration.

Polish Mode erases surface stains and enhances natural tooth whiteness.

Massage Mode provides gentle gum stimulation during brushing to soothe sensitivity.

9. After selecting your mode, set the intensity level by pressing the function selection button. The MORI 2.0 offers 5 intensity levels (10, 30, 50, 70, 90) for every mode. Note that the brush features splash prevention technology: when activated, it initiates low-intensity vibrations to prevent toothpaste splatter. Full-power operation begins only when the pressure sensor detects contact with teeth.

10. The MORI 2.0 features a memory function: after 10 seconds of use in any mode, it saves the setting and automatically reactivates it during your next brushing session.

11. To remind you to switch brushing quadrants, the MORI 2.0 pauses briefly during operation. For a 2-minute brushing session, it pauses every 30 seconds; for 3 minutes, every 45 seconds; for 4 minutes, every 60 seconds. If "Unlimited" brushing time is selected, no pauses will occur.

12. The MORI 2.0 sonic toothbrush features a pressure sensor. A visual indicator on the rear panel monitors brushing pressure through color-coded alerts, vibration patterns, and an audible "beep-beep-beep" signal.

Blue light: Insufficient pressure

Green light: Optimal pressure

Red light: Excessive pressure

13. The toothbrush will automatically power off after completing the preset brushing time. To manually stop brushing earlier, briefly press the power button after 5 seconds of operation.

14. After brushing completes, the display will show both the remaining battery level and your cleaning results with a 3D oral map.

Cleaning Effectiveness Color Codes on the Map:


White: Area fully cleaned — brushing duration met recommended standards.

Blue: Insufficient cleaning time — this oral zone was partially cleaned.


Orange: Uncleaned area — brushing wasn't initiated here.

15. After 3 months of use (based on twice-daily brushing: 2 uses/day × 30 days/ month × 3 months = 180 sessions), a brush head replacement reminder will appear on the display when you power on the toothbrush.

Usage Method

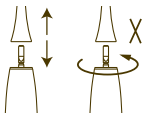1. Install the brush head.

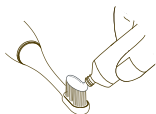
2. Wet the brush head with water and apply the recommended amount of toothpaste.


3. Place the brush head in your mouth. Position it at a 45-degree angle toward the gums. Gently place the bristles at or slightly above the gum line.

Note: The middle section of the brush head must maintain constant contact with your teeth.

4. Turn on the toothbrush and select a convenient cleaning mode.
5. Begin brushing using gentle sweeping motions with the brush head, moving from the gumline toward the chewing edge of each tooth.

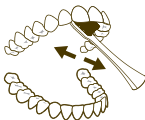
6. For cleaning the inner surfaces of upper teeth, position the toothbrush vertically. This aligns the brush head perpendicular to each tooth for optimal cleaning.

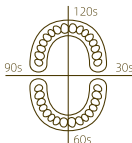7. Dentists recommend brushing for at least 2 minutes. Divide your mouth into 4 quadrants, spending 30 seconds cleaning each section. The MORI 2.0 offers customizable brushing durations (2/ 3/ 4 minutes or Unlimited), with quadrant timing automatically adjusted to your selected duration. To support this technique, the brush pauses briefly between quadrants — reminding you to move to the next area.

8. Review your brushing results. The MORI 2.0 divides all teeth into 16 zones — including anterior, palatal/ lingual, and occlusal surfaces — displayed in a dental map diagram. During brushing, the device identifies cleaned and missed areas, marking them on the 3D map to help you achieve comprehensive cleaning coverage.

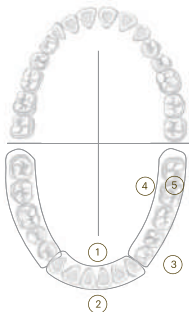1. Inner side of incisors
2. Front side of incisors
3. Outer part of large molar teeth
4. Inner part of large molar teeth
5. Occlusal (bite) surface of large molar teeth

Cleaning

- Brush heads of the MORI 2.0 electric sonic toothbrush: Rinse thoroughly under running water after every use and allow the brush head to air dry naturally.

- Cleaning the Toothbrush Body

The brush body may be cleaned directly with water; afterward, allow the device to air dry naturally. Do not use water or acidic cleaning agents at temperatures exceeding 40°C.

The MORI 2.0 kit includes a UV sanitization case for brush heads.


Technical Specifications

Parameter	Value
Power	2 W
Operating Voltage	5 V
Power Source	Lithium-ion Battery
Battery Capacity	600 mAh
Sterilization Efficiency	>99.3%
UV Wavelength	275 nm
Case Capacity	2 brush heads
Smart Timer	5 minutes
Operating Time	200 min
Charging Time	3 hours

Purpose


Residual toothpaste and saliva accumulate on bristles after use, creating a moist environment where bacteria rapidly multiply. The UV sterilizer included with the MORI 2.0 toothbrush eliminates bacteria from brush heads and prevents future microbial buildup. This provides a quick, efficient disinfection solution.

Operation and Usage Method

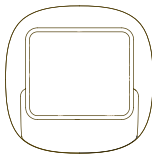
1. The sanitization case is designed for wall mounting, equipped with a holder featuring adhesive backing for installation on any vertical surface. To mount it: pull the holder downward to detach it from the case body, remove the protective film from the adhesive surface, affix the holder to your wall or other vertical surface, then insert the case into the holder's grooves and pull downward to secure it.

2. Insert the brush heads into the case: position the bristle head inside the closed compartment while allowing the handle to extend beyond the case.

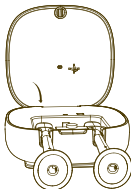3. Close the case lid to initiate operation.

4. Sterilization begins automatically, completes after 5 minutes, and the UV lamp switches off independently.

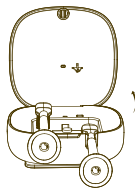5. Open the case and remove the brush heads.

Cleaning

After use, open the sanitizer lid to air out the compartment. Regularly wipe the case exterior with a damp cloth.

Storage

Store the device in a well-ventilated area away from direct sunlight.

Keep out of reach of children.

Charging

1. When the battery is low, the blue indicator light on the case blinks.
2. Insert the Type-C end of the charging cable into the case and the USB end into a power adapter to begin charging.
3. The indicator will glow steadily when fully charged.

Precautions

1. Do not attempt to repair the device yourself.
2. Avoid using the unit in high-temperature environments.
3. Never place foreign objects, explosive materials, or flammable substances inside the case.
4. Do not open the device while operational.

Storage & Safety

- Avoid dropping

Stop if the toothbrush body or any brush head appears damaged.

- Store away from heat sources.
- Do not disassemble the MORI 2.0 sonic toothbrush or UV sanitization case.

For repairs, contact the manufacturer-authorized service center listed in the warranty section.

Charging

1. When the battery level falls below 20%, a low charge notification will appear on the toothbrush display after powering on. Press the function selection button to dismiss this alert.


2. At battery levels below 10%, the display will turn off and the toothbrush will automatically shut down 5 seconds after being powered on.

3. Before charging, ensure the brush body is completely dry. Wipe it thoroughly if water residue is present.

Important: Avoid charging in high-humidity environments and prevent water contact with the brush body or charging cable during charging!

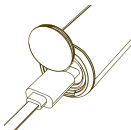4. Insert the USB end of the cable into a power adapter, and the Type-C end into the recharging port located at the base of the brush body.

5. If the cable is properly connected, the toothbrush display will light up.

6. During charging, the battery level will be displayed in real time on the screen.

7. Wait until 100% charge is reached, then disconnect the toothbrush from the power adapter.

Precautions

The MORI 2.0 toothbrush may be used by children over 8 years of age and persons with sensory, physical, or mental limitations, provided they have received instruction on the proper use of this device and are aware of safe operation methods and potential hazards of misuse.

It is prohibited to use this device as a toy.

Brushing and device usage by children must occur only under adult supervision.

Disposal

The MORI 2.0 is manufactured and assembled using high-quality materials and components suitable for recycling and reuse. When disposing of this product, please do not discard it with household waste. Care for the environment by recycling this electrical appliance according to local regulatory requirements.

Warranty

A 1-year warranty applies to the WAU MORI 2.0 Electric Sonic Toothbrush. During the warranty period, the product shall be repaired or replaced if manufacturing defects are confirmed. This warranty does not cover malfunctions caused by improper use or transportation, nor products subjected to unauthorized repair.

MORI 2.0

ELECTRIC SONIC TOOTHBRUSH

Model: **RM-T40**

Material: **Polycarbonate plastic**

Dimensions (H×W×D): **180 × 26 × 26 mm**

Weight: **110 g**

Voltage: **~3.7 V**

Power: **5 W**

Made in: **PRC**

Authorized manufacturer representative: I. P. Umantsev A. N.
Leningradsky Prospekt 37, Office 14-14, Moscow 125167, Russia
info@wauglobal.com

Important Notice: To ensure continuous product improvement, the manufacturer reserves the right to modify design and technical specifications. The latest version of this manual is available at: wauglobal.com